





imp De Mayo Knee Positioner® Instructions For Use

1. Replace OR Table Pads with green *imp* pads.


Note: It is not necessary to remove the X-Ray cassette plates.

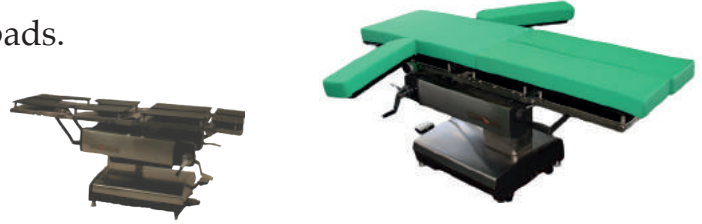


2. Position patient with gluteal fold at separation of pads.

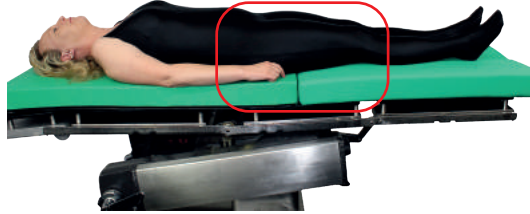

3. During prep of surgical leg, remove single OR Table Pad.

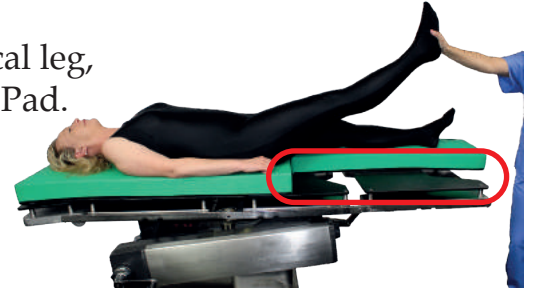

4. Create a 3" - 5" fold in the final drape under the buttocks.

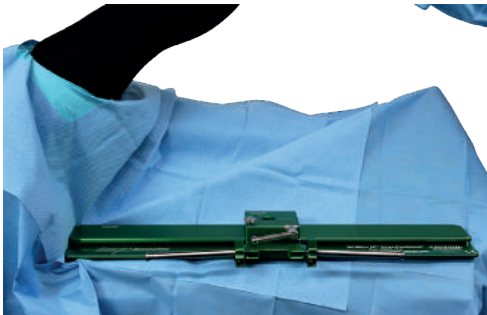

5. Place De Mayo Knee Positioner® in the well to achieve maximum flexion.




6. Check flexion for final approval of freedom of the drape to allow full range of motion when the Positioner is locked on drapes.



imp® De Mayo Knee Positioner® Instructions For Use

7. Position clamp directly under gliding pins and lock to the OR Table Side Rail utilizing either

De Mayo Single Lever Clamp ~ or ~ De Mayo Push Button™ Clamp ~ or ~ De Mayo Push Button™ Clamp-Pin Locking

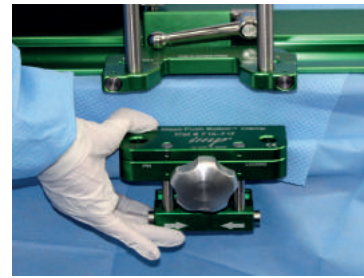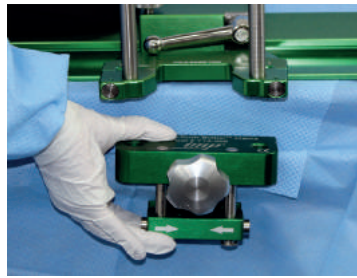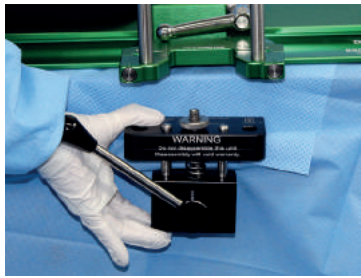

PLEASE REFER TO APPLICABLE CLAMP INSTRUCTIONS FOR USE

8. Protect your patient from potential Pressure Ulcers. Insert the *imp* Patient Protective Pad into the sterile boot.




STERILE EO



9. Wrap cohesive bandage around the foot starting with a minimum of six (6) foot wraps, tear, and finish the wrap above the bracket around the calf.



Avoid wrapping cohesive bandage over the bracket on the back of the boot.

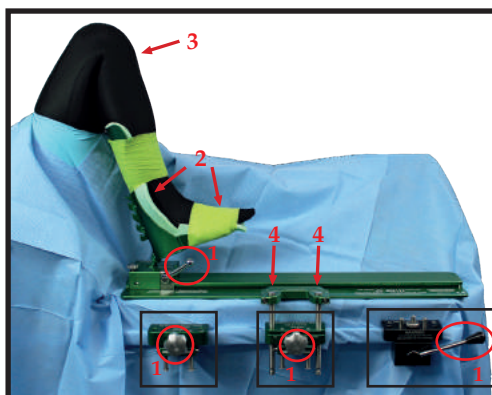

Set Up Test Procedure



TIGHTEN CLOCKWISE TO LOCK KNOBS AND LEVERS:

- 1 ~ Check all knobs and levers are locked
- 2 ~ Patient is fully protected by imp® foam and cohesive wrap
- 3 ~ Patient's leg can reach full flexion
- 4 ~ Check that both sliding pins are fully seated

imp® Innovative
Medical Products, Inc.

87 Spring Lane, PO Box 8028
Plainville, CT 06062

PH: (860) 793-0391 FAX: (866) 459-1805

Toll Free PH: 1-800-467-4944

www.impmedical.com

